

DATASHEET

Lenti-Mkate2-P virus

Cat. VSL-0049

Titer:	~10 ⁶ IU/ml
Medium:	DMEM, 10% FBS
Volume:	10 ml (10x1ml)
Selection marker:	Puromycin
Storage:	-70°C
How to use:	Thaw the recombinant lentivirus supernatant in a 37°C water bath; remove it from the bath immediately when thawed.

Description: Ready-to-use virus expressing far-red fluorescent protein mKate2 can be used as a control in experiments with lentivirus transductions of target cells or for the generation of stable cell lines.

mKate2 is the next generation of monomeric far-red fluorescent protein TagFP635 (mKate) [Shcherbo et al., 2007; Shcherbo et al., 2009]. Lenti-Mkate2 virus is VSV-G pseudotyped lentivirus that can be used as a control for infection of target cells. It was generated by co-transfection of 293T packaging cells with pLVmKate2 plasmid (Fig.1) and Packing Mix. Virus-containing supernatant has been harvested 48-72 hours post-transfection, centrifuged at 3000rpm for 15 minutes at +4°C to pellet debris, and then filtered through 0.45µm PVDF syringe.

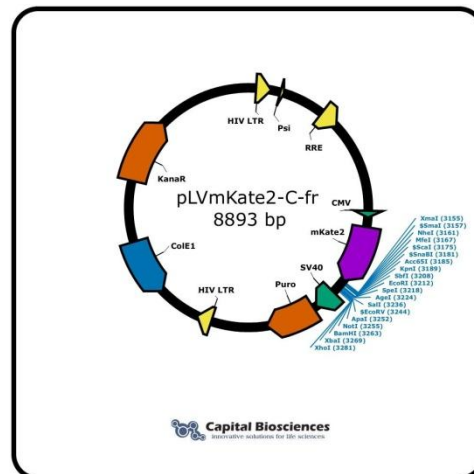


Fig.1. Map of pLVmKate2 vector

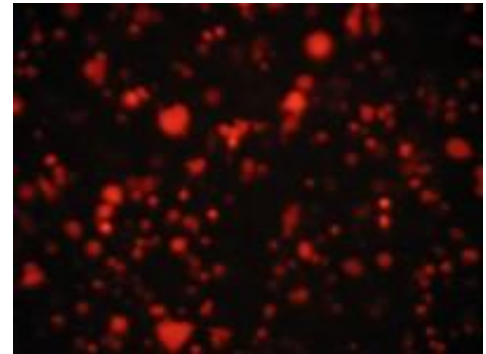


Fig.2. 293 cells infected with Lenti-mKate2 virus

Notice to Purchaser: All products containing mKate2 are intended for research use only. The Products are covered by U.S. Pat. # 7,417,131 and other Evrogen Patents and/or Patent applications pending. By use of these Products, you accept the terms and conditions of the applicable Limited Use Label License.

For research use only

9700 Great Seneca Hwy, Suite 247
Rockville, MD 20850
Tel: 800-475-2812; Fax: 240-363-0739